COMPRESSOR NEBULIZER

CND-501S

SAMSUNG AMERICA, INC.

105 CHALLENGER ROAD 6th FLOOR
RIDGEFIELD PARK, NJ 07660

Printed in China

Toll Free Customer Service

1.800.261.0130

www.SAMsungHealthyLiving.com

SAMSUNG

SAMSUNG AMERICA

Your Healthy Living® Compressor Nebulizer by Samsung America is intended for use in the treatment of asthma, COPD and other respiratory ailments in which aerosolized medication is required during therapy. Please consult with your physician and/or pharmacist to determine if your prescription medication is approved for use with this Nebulizer. This Nebulizer is intended for single patient use.

Federal United States law restricts this device to sale by, or on the order of a licensed physician.

**Instruction
Manual**

Available in English and Spanish

Contents

Please read instruction manual carefully before operating.

Introduction	2
Product Identification	2
Important Safeguards	3
Operation of the Unit	4
Assembling your Nebulizer	
Operating your Nebulizer	
Important Tips	
After Every Use	
Cleaning and Maintenance	6
Cleaning Procedures	
Air Filter Replacement	
Troubleshooting	7
Specifications	8
Accessories	
Warranty	9
Spanish Version	11

For your future reference -

Date of Purchase: _____

Serial Number: _____

Attach your receipt here:

Introduction


Thank you for choosing the Healthy Living® Compressor Nebulizer by Samsung America. With proper care and use, it will provide you with many years of reliable treatment.

This product is developed for the successful treatment of asthma, COPD, allergies and other respiratory ailments in which aerosolized medication is required. It creates a stream of air that travels through a clear tube to the nebulizer. When air enters the nebulizer, it will convert the prescribed medication into aerosol mist for easy inhalation. Treatments are delivered quickly, safely and conveniently making this unit ideal for all ages.

This Nebulizer should be used under the supervision of a licensed physician and/or a respiratory therapist. We encourage you to thoroughly read this instruction manual to learn about the features of this product. Any use of this product other than its intended use should always be avoided.

Product Identification

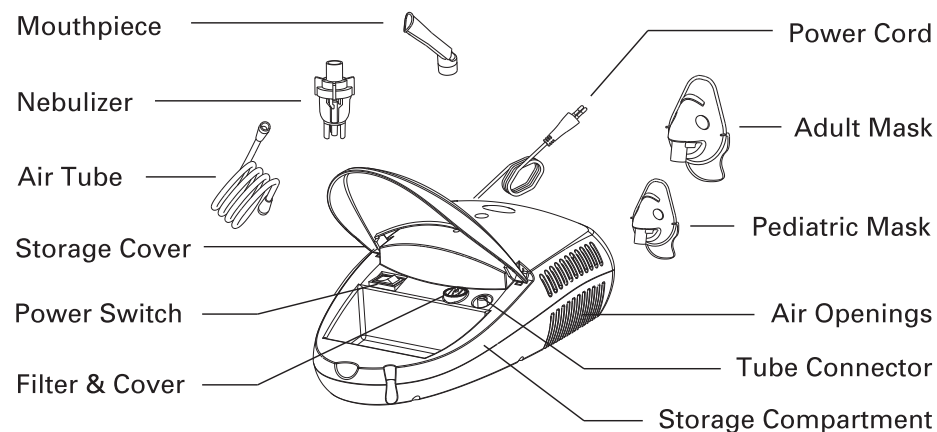Contents: Instructional DVD, Storage Case, Warranty Card

Important Safeguards

READ ALL INSTRUCTIONS BEFORE USING

Failure to read and observe all precautions could result in personal injury or equipment damage.

PRECAUTIONS:

To avoid electrical shock:

- Keep the unit away from water.
- Do not immerse the power cord or the unit into water or other liquid.
- Do not use while bathing.
- Do not reach for a unit that has fallen into water – immediately unplug the unit.
- Always unplug the unit immediately after using.

OPERATING CAUTIONS:

- Adult supervision is highly recommended when unit is used by children.
- Never operate the unit if it has any damaged parts (including power cord or plug), if it has been dropped or submersed in water. Promptly send it to a service center for examination and repair.
- The unit should not be used where flammable gas, oxygen or aerosol spray products are being used.
- Do not tilt or shake the unit when in use.
- Disconnect the unit from the electrical outlet before cleaning, filling and after each use.
- Do not use other attachments unless recommended by manufacturer.

STORAGE CAUTIONS:

- Do not store the unit under direct sunlight, high temperature or humidity.
- Keep the unit out of reach of small children.
- Always keep unit unplugged when not in use.

CLEANING CAUTIONS:

- Before each use, check the air filter, Nebulizer, mouthpiece, and any other optional components.
- Dirty or worn parts should be replaced.
- Do not immerse the unit in water. It may damage the unit.
- Before cleaning, disconnect the unit from the electrical outlet.
- After each use, clean all necessary parts as instructed in this manual.

TO REDUCE THE RISK OF PERSONAL INJURY:

1. Connect this unit to an appropriate voltage outlet for your model.
2. Do not run the unit unattended.
3. If any abnormality occurs, discontinue use immediately until the unit has been examined and repaired.
4. Never block the air openings of the main unit or place it where the air openings may be obstructed.
5. Always consult with your physician to determine if your prescription medication is approved for use with this compressor nebulizer.

PLEASE SAVE INSTRUCTION MANUAL FOR FUTURE REFERENCE

Operation of the Unit




Follow the cleaning instructions in this manual under “Cleaning & Maintenance” prior to using your nebulizer for the first time or after it has been stored for an extended period of time.

Remember: Always unplug the Compressor and make sure the POWER switch is turned to the “off” position before cleaning, assembling and before or after each use.

Assembling your Nebulizer


- 1) Place the Compressor on a flat, stable surface within reach, with the handle facing towards you.
- 2) Open the storage compartment by lifting the lid. Remove the contents (Nebulizer, mouthpiece and air tube).
- 3) Gently twist and pull straight up on the lid of the Nebulizer to separate into two parts (medicine cup and cover).
- 4) Make sure that the purple jet is properly installed inside the medicine cup. (The stem inside the medicine cup inserts into the tube of the purple jet)
- 5) Add the prescribed amount of medicine to the medicine cup.
- 6) Reassemble the Nebulizer by carefully twisting the medicine cup and cover together. Make sure that the two parts fit securely.

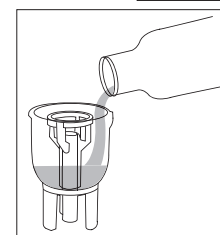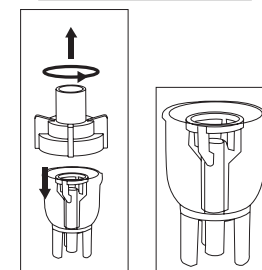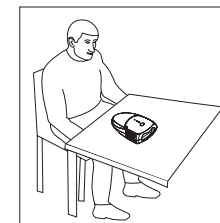Operating your Nebulizer

Note: The Nebulizer is operable at up to a 45° angle.

- 1) Attach one end of the air tube connector to the air stem.

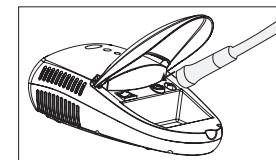
2) Carefully attach the opposite end of the air tube connector to the stem at the base of the Nebulizer.

3) Attach the angled mouthpiece or mask to the top of the Nebulizer.

If using an optional tee adapter set:

- a) Attach the tee adaptor to the top of the Nebulizer.
- b) Connect one end of the tee adapter to the straight mouthpiece and the other to the extension tube.

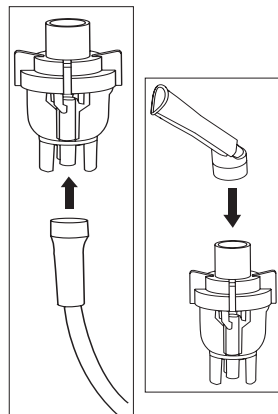
Note: The tee adapter is not needed when using a mask.

4) Check that the POWER switch is in the “off” position. Plug the power cord into a main outlet.

5) Press the POWER switch to the “on” position.

6) Start your treatment according to your physician’s or therapist’s instructions.

Note: The crescent-shaped groove on the side of the compressor can be used to support the nebulizer during treatments.

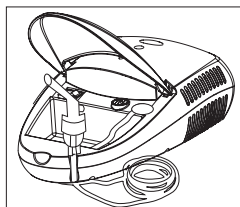7) Once the treatment is complete, press the POWER switch to the “off” position and unplug the power cord from the outlet.

IMPORTANT TIPS!

The Compressor will automatically shut off if it becomes overheated. If/when this happens, immediately:

- 1) Press the POWER switch to the “off” position.
- 2) Unplug the power cord from the outlet.
- 3) Allow the motor to cool for 30 minutes.

Before restarting the unit, make sure that the air vents are not obstructed.

After Every Use

- 1) Unplug the unit from the power source.
- 2) Allow the unit to completely cool.
- 3) Carefully detach the air tubing from the Nebulizer and pour out any remaining medication.
- 4) Follow the cleaning procedures provided in this manual.

Cleaning and Maintenance


Cleaning Procedures

Please consult with your physician or therapist prior to cleaning regarding the use of commercial disinfectants and appropriate cleaning methods.

Cleaning Technique

(Performed after each treatment):

- 1) Make sure that the POWER switch has been turned to the “off” position and the unit has been disconnected from the power source.
- 2) Disconnect the air tube from the base of the Nebulizer.
- 3) Gently twist and pull up on the cover of the Nebulizer to open and separate.
- 4) Wash all parts of the Nebulizer, the mouthpiece, and other optional components in warm water with a mild detergent.
- 5) Rinse the Nebulizer and components with hot tap water.
- 6) Allow all parts to completely air dry on a clean paper towel.

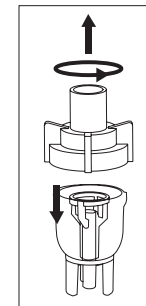Note: The air tube does not have to be washed as only filtered air passes through it. Keep the outer surface of the air tube dust-free by wiping regularly.

Disinfecting Technique

(Performed after the last treatment of the day or before first time use):

- 1) Follow steps 1-3 under “Cleaning Technique”.
- 2) Submerge the Nebulizer parts in a clean container or bowl that contains a vinegar/distilled water solution of one part white vinegar to three parts distilled water. Let the parts soak for 30 minutes.
- 3) Follow steps 5-6 under “Cleaning Technique”. Do not towel dry Nebulizer parts or reuse disinfecting solution as this could cause contamination.

Note: For the Nebulizer only, clean in the dishwasher using the top shelf. Do not place any Nebulizer parts in the dishwasher. Doing so may cause damage.

Cleaning the Compressor


Wipe the compressor daily using a soft, slightly water-moistened cloth.

Note: Do not place Compressor in dishwasher. Any other form of cleaning or cleaning agents may damage the finish on the unit.

Air Filter Replacement

It is important to change the air filter approximately every 30 days or when the air filter turns grey.

- 1) Remove the air filter cover by gently pulling forward.
- 2) Discard the gray air filter.
- 3) Replace with a new clean air filter.
- 4) Securely reattach the air filter cover to the unit.

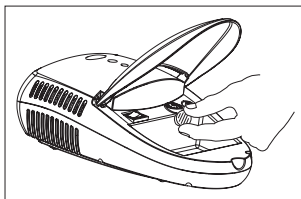Note: Air filters cannot be cleaned or washed. Do not use cotton or any other materials. Only use Healthy Living® filters. Do not operate without a filter.

Troubleshooting

If any abnormality occurs during use, please check and correct the following:

Condition	Correction
Unit does not operate when POWER switch is pressed.	Check the AC connection to the outlet
No misting or low rate of misting	Check that there is medication in the Nebulizer cup
	Check the Nebulizer for damage
	Check the position of the purple jet inside of Nebulizer
	Make sure that air tube and other components are properly attached
	Check the air filter and replace if necessary
	Check the air tube and replace if necessary

Specifications

**HEALTHY
LIVING®**

Power	120V, 60Hz
Power Consumption	80VA
Medication Capacity	5 ml
Particle Size	0.5 to 10um
MMAD	4um
Sound Level	< 60 dBA
Maximum Capacity of Medication	13 ml (cc)
Average Nebulization Rate	Min 0.2 ml/min
Compressor Pressure Range	30 to 36 psi
Operating Pressure Range	8 to 16 psi
Liter Flow Range	6 ~ 8 lpm
Operating Temperature Range	50°F to 104°F (10°C to 40°C)
Operating Humidity Range	10 to 95% RH
Storage Temperature Range	-4°F to 158°F (-20°C to 70°C)
Storage Humidity Range	10 to 95% RH
Dimension (L x W x H)	11.4" x 7.2" x 4.2" (290 x 182x 106mm)
Weight	4.4 lbs (2.0 kg)
Standard Accessories	Nebulizer, Air Tube, Adult & Pediatric Masks, Mouthpiece, Filters (5 pcs)

For purposes of improvement, specifications are subject to change without notice.

Protection against electric shock:

- ☐ Class II Equipment
- ☒ Type BF applied parts: Mouthpiece & Masks

Protection against harmful ingress of water: IPX0 (Ordinary)

Degree of safety in the presence of flammable anesthetics or oxygen: No AP/ APG (Not suitable for use in the presence of flammable anesthetics or oxygen)

Mode of operation: continuous

With respect to electric shock, fire and mechanical hazards only in accordance with UL 1431

Accessories

CNK-510S (Replacement Nebulizer Kit): Nebulizer, and Air Tube

CNF-520S (Replacement Filters): 10 pcs

Warranty

Healthy Living® Compressor Nebulizer Limited 5 Year Warranty

Samsung America, Inc. warrants to the original purchaser of the nebulizer that your Healthy Living® Compressor Nebulizer will be free from defects in material and workmanship for 5 years from the date of purchase, under normal use. This product is for home use only and is not intended for commercial or clinical use. If during this 5 year period, the compressor nebulizer does not work properly because of defect in materials or workmanship, Samsung America, Inc. will replace it with a new Healthy Living® Compressor Nebulizer or equivalent product free of charge (except you must pay the shipping charges set forth below). Purchase receipt or other proof of date of original purchase is required before warranty performance.

This warranty extends solely to failures due to defects in materials or workmanship occurring during normal use. It does not cover normal wear of the product. This warranty is void if: the product housing has been removed, if the product's label, logo or serial number have been removed or if the product fails to function properly as a result of accident, misuse, abuse, neglect, mishandling, misapplication, faulty installation, setup, adjustments, improper maintenance, alteration, maladjustment of controls, modification, power surges, commercial or clinical use of the product, service by anyone other than an authorized service center or acts beyond the manufacturer's or distributor's control.

The warranty and remedy provided herein are exclusive and in lieu of all other express and implied warranties and unless stated herein, any statements or representations made by any other person or entity are void. NEITHER THE MANUFACTURER, DISTRIBUTOR NOR THEIR AFFILIATES SHALL BE LIABLE FOR ANY LOSS, INCONVENIENCE, INJURY, OR DAMAGE INCLUDING DIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE, IMPROPER USE, FAILURE TO HEED ANY WARNING, DISCLAIMER OR INSTRUCTION OR INABILITY TO USE THIS PRODUCT, WHETHER ARISING OUT OF ANY BREACH OF THIS WARRANTY OR ANY OTHER LEGAL THEORY ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING THE WARRANTIES OF MERCHANTABILITY, ARE LIMITED TO THE APPLICABLE WARRANTY PERIOD SET FORTH ABOVE. Manufacturer and Distributor assume no responsibility for injuries, damages or penalties incurred or suffered resulting from the use of this product in an improper manner or manner or location other than for that for which it is intended or from failure to heed any warning, disclaimer or instruction.

Some states do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so the above exclusions or limitations may not apply to you. This warranty gives your specific legal rights and you may also have other rights which vary from state to state.

Should your product prove defective during this warranty period, please return the product in the original carton prepaid to:

SAMSUNG AMERICA, INC.
c/o Healthcare Division
105 Challenger Road, 6th Floor
Ridgefield Park, NJ 07660
Attn: Customer Service

Include **\$8.95** for return shipping and insurance. (CHECK OR MONEY ORDER ONLY PLEASE, NON-REFUNDABLE).

Include with your nebulizer a letter stating: your name, street address, daytime phone number, detailed description of the problem, a copy of your dated cash register receipt or other proof of purchase.

Sorry, we cannot deliver to P.O. Boxes or outside the U.S.

This warranty applies only to products purchased in the U.S. For products purchased in, but used outside, the U.S., this warranty covers only warranty services within the U.S. (and does not include shipping outside the U.S.).

© 2005 Samsung America, Inc. All rights reserved.

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